


Capillary HPLC:


PLRP-S Separation of Parabens (Preservatives)

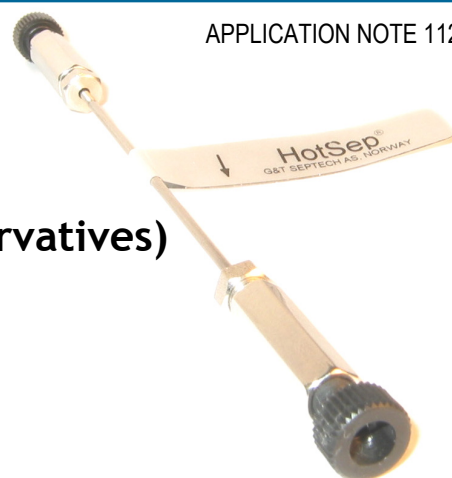

App No 106

Parabens

Column: 3 μ PLRP-S 100Å
Dimensions: 0.3 x 150 mm
Order No: S-161-0315
Mobile Phase: MeCN-water (40:60)
Flow Rate: 5 μ L/min
Detector: UV @ 254 nm

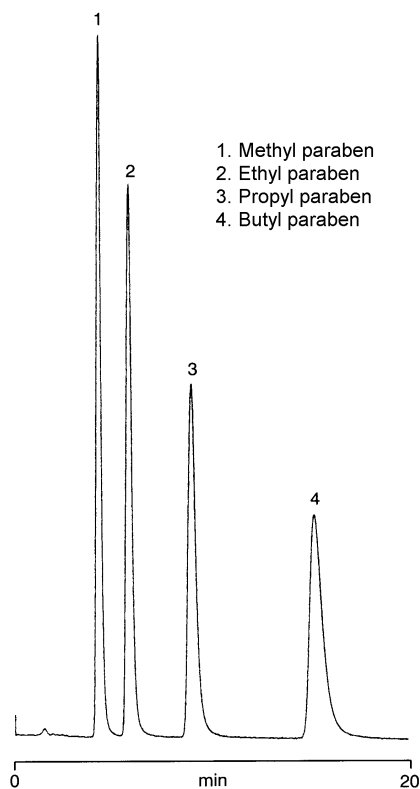

Comments:

Parabens are a group of chemicals widely used as preservatives in food and cosmetic and therapeutic products. Some parabens are also found at low levels in nature. Paraben is the common name for this class of chemical, however, they are also known by other names such as esters of p-hydroxybenzoic acid.

ORDERING INFORMATION: HotSep® PLRP-S, 3 μ , 100Å

ID	3 cm	5 cm	10 cm	15 cm	Guard	Guard/5pk
1.0 mm	S-161-1003	S-161-1005	S-161-1010	S-161-1015	G-161-10-1	G-161-10-5
0.5 mm	S-161-0503	S-161-0505	S-161-0510	S-161-0515	G-161-05-1	G-161-05-5
0.3 mm	S-161-0303	S-161-0305	S-161-0310	S-161-0315	G-161-03-1	G-161-03-5
0.1 mm	---	---	S-161-0110	S-161-0115	---	---
75 μ m	---	---	S-161-00710	S-161-00715	---	---