

## Microbore HPLC:


### Fast separation of fluorene (PAH) on 50 mm long column




App No 133

#### Fluorene


**Column:** 5µ Kromasil Cyano 60Å  
**Dimensions:** 1 x 50 mm  
**Order No.:** S-124-1005  
**Mobile Phase:** MeCN-water (70:30)  
**Flow Rate:** 35 µL/min  
**Temp.:** 25C  
**Detector:** UV @ 254 nm

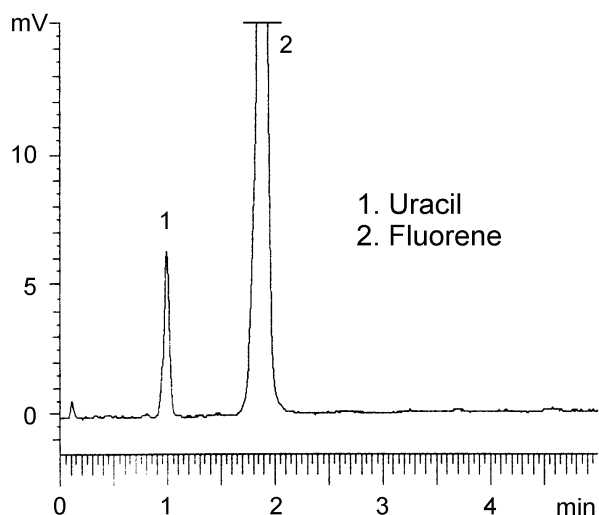

1.

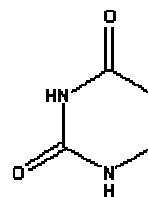

2.

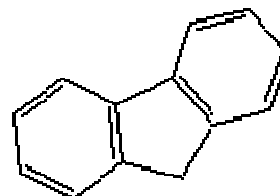

#### ORDERING INFORMATION: HotSep® Kromasil Cyano (CN), 5µ, 100Å

ID	5 cm	10 cm	15 cm	25 cm	Guard	Guard/5pk
1.0 mm	S-124-1005	S-124-1010	S-124-1015	S-124-1025	G-124-10-1	G-124-10-5
0.5 mm	S-124-0505	S-124-0510	S-124-0515	S-124-0525	G-124-05-1	G-124-05-5
0.3 mm	S-124-0305	S-124-0310	S-124-0315	S-124-0325	G-124-03-1	G-124-03-5
0.1 mm	---	---	---	---	---	---
75 µm	---	---	---	---	---	---